

T-Shaped and L-Shaped Spreaders




T-Shaped Plastic Spreaders

The Sterile T-Shaped Spreaders are designed for spreading and dispersing liquids across the surface of agar culture plates. Spreaders are used for performing bacterial CFU plate counts on water, milk, urine and other liquid samples. Copan T-Shaped Spreaders have all smooth rounded surfaces to prevent any cutting or gouging of agar during inoculation of culture plates. The T-Shape design means that even pressure is applied across the entire length of the T-bar in contact with the agar surface. T-Shaped Spreaders are available individually sterile wrapped or in packs of 5 and 10 pieces in a peel pouch.

L-Shaped Plastic Spreaders

for spreading liquid samples The foot of the spreader has a completely smooth rounded surface, free of rough edges and imperfections in the plastic. This feature enables even spreading of liquid samples across the surface of agar plates without gouging or cutting the medium. L-Shaped Spreaders are available in two sterile pouch sizes - single wrapped or five spreaders per pack.

Catalog No.	Description	Packaging	Case
173CS01	T-Shaped Spreader	1 / peel pack	500 PCS
173CS05	T-Shaped Spreader	5 / peel pack	1000 PCS
173CS10	T-Shaped Spreader	10 / peel pack	1000 PCS
174CS01	L-Shaped Spreader	1 / peel pack	500 PCS
174CS05	L-Shaped Spreader	5 / peel pack	1000 PCS

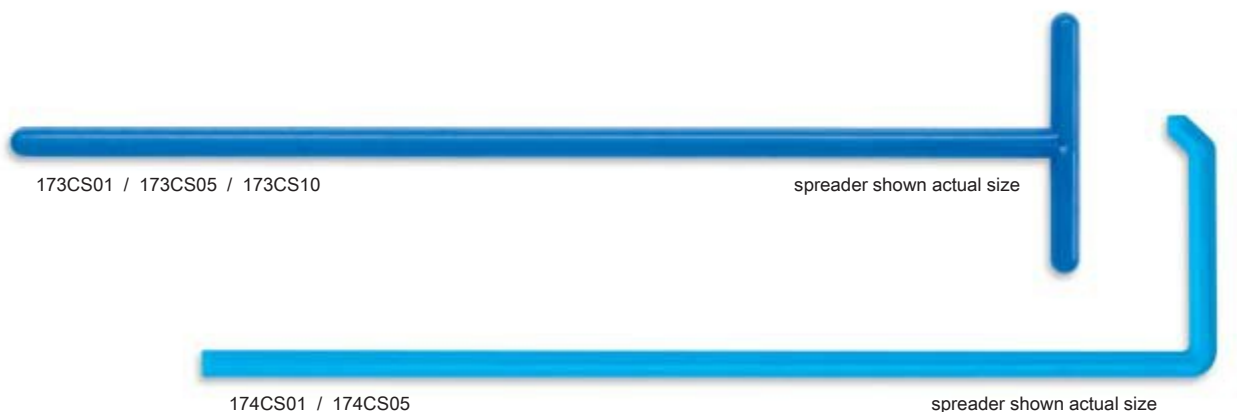

BioGenerica
Your Partner for Life Science and Research Products

BIOGENERICA S.R.L.
Via Pescara 7 Int. 23 - 95030 Mascali (CT) - Italia
Tel. +39 095914678 - Fax +39 0952246205 - info@biogenerica.it
www.biogenerica.it